

ceveron®
MFU 500



MICROPARTICLE FILTRATION UNIT

The new Ceveron® MFU 500 enables easy, quick and standardised removal of MICROPARTICLES from platelet poor plasma.

EASY & QUICK



Pipette sample into the filter plate



Push the "START" button and filtration starts automatically



After filtration is terminated push the "VENTILATION" button



Filtered sample can be collected from the receiver plate

SMALL & COMPACT




Filtration is as effective as ultra-centrifugation to remove thrombogenic MP's from PPP.

Dimensions MFU 500:

Length: 26.00 cm
Width: 20.00 cm
Height: 16.05 cm
Weight: 3.58 kg

STANDARDIZED

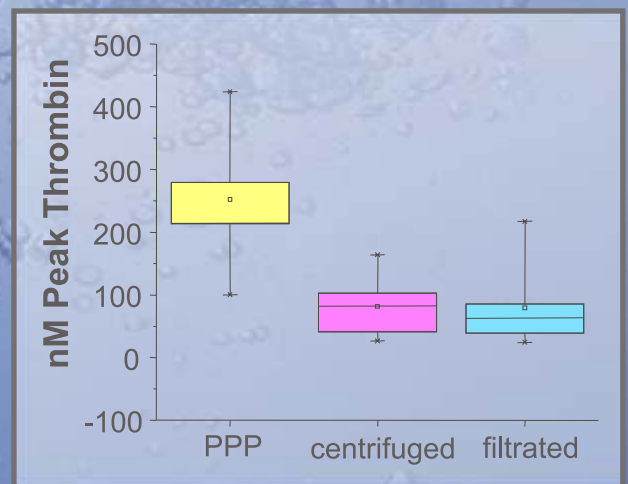

Technoclone GmbH
Brunner Str. 59/5
1230 Vienna, Austria
Tel: +43 1 86373-10
Fax: +431 86373-44

Technoclone Deutschland GmbH
Hans Bunte Str. 8-10
69123 Heidelberg, Germany
Tel: +49 6221 830 441
Fax: +49 6221 830 437

Technoclone Ltd.
Eclipse House 7, Curtis Road,
Dorking/Surrey RH4 1 EJ
Tel: +44 1306 888 777
Fax: +44 1306 883 883



technoclone

...our quality is your success!