

TECHNOCLOT PT-PLUS

PRODUCT DESCRIPTION

- ◆ Technoclot PT Plus is a standardized thromboplastin obtained of rabbit brain
- ◆ ISI-value of about 1.1
- ◆ Each batch of Technoclot PT Plus is standardized according to ISTH guidelines using the International Reference Preparation as reference
- ◆ Technoclot PT Plus is insensitive to heparin concentrations up to 2 I.U. / mL in plasma
- ◆ Optimized reaction curves allow test performance on all mechanical and optical based coagulation analyzers with very low CVs
- ◆ High sensitivity to FII, FV, FVII and FX




APPLICATION





Technoclot PT Plus is indicated

- ◆ as screening test for bleeding tendencies due to the extrinsic coagulation system
- ◆ as liver function test
- ◆ for determination of coagulation factors II, V, VII and X
- ◆ monitoring and adjustment of oral anticoagulant therapy



HEPARIN INSENSITIVITY

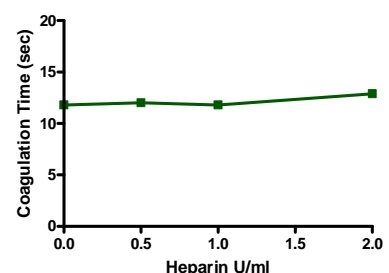

FACTOR SENSITIVITY

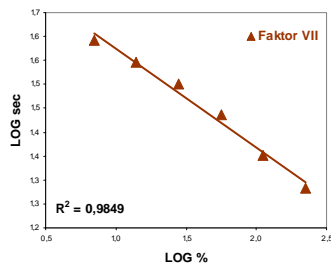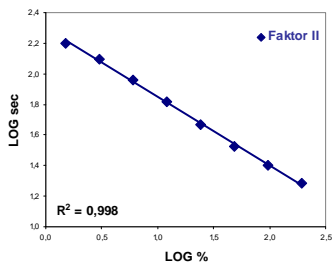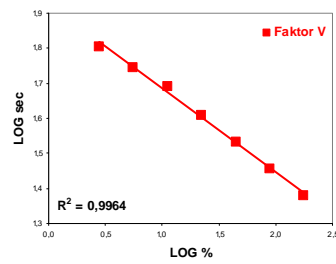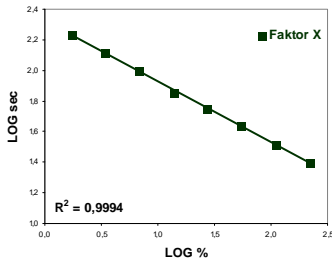

GOOD CORRELATION WITH OTHER PT REAGENTS

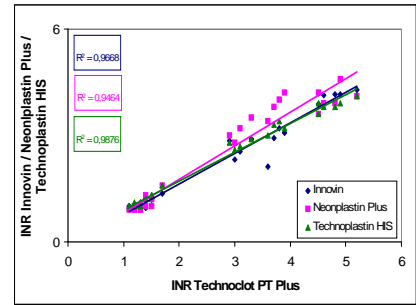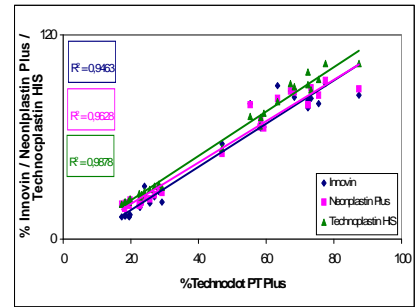

ORDER INFORMATION

REF 5003220 Technoclot PT PLUS 5 x 4 mL

REF 5003222 Technoclot PT PLUS 10 x 10 mL

Calibrators

REF 5010004 AK Calibrant 4 x 1 mL

Controls

REF 5011050 Coagulation Control AK 5 x 1 mL

REF 5011060 Coagulation Control AK 50 x 1 mL

REF 5020040 Coagulation Control N 5 x 1 mL

REF 5020050 Coagulation Control N 50 x 1 mL

REF 5021055 Coagulation Control A 5 x 1 mL

REF 5021060 Coagulation Control A 50 x 1 mL

REF 5021070 Technoclot Control A 10 x 1 mL

REF 5021075 Technoclot Control A 50 x 1 mL

REF 5020040 Technoclot Control N 5 x 1 mL

REF 5020040 Technoclot Control N 50 x 1 mL