

technoclone

# CEVERON<sup>®</sup> MFU 500

Microparticle filtration unit






*A clear view into the future!*

**Ceveron<sup>®</sup> MFU 500**  
 enables easy, quick and standardised  
 removal of *MICROPARTICLES* from  
 platelet poor plasma.



**EASY & QUICK**

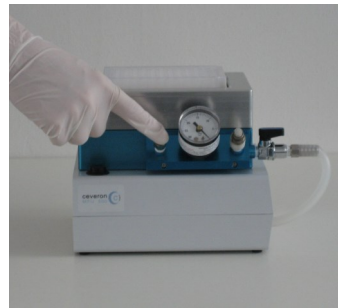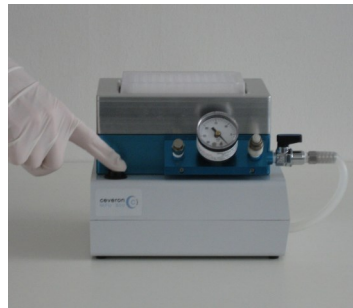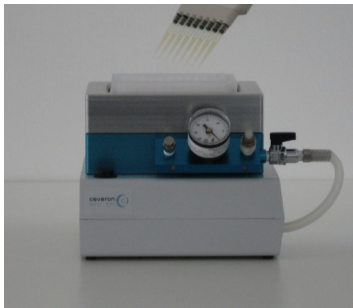

Pipette sample into  
 the filter plate

Push the "START"  
 button and filtration  
 starts automatically

After filtration is  
 terminated push the  
 "VENTILATION" button

Filtered sample can be  
 collected from the  
 receiver plate


**SMALL & COMPACT**



Filtration is as effective as  
 ultra-centrifugation to re-  
 move thrombogenic MP's  
 from PPP.

**Dimensions MFU 500:**  
 Length: 26.00 cm  
 Width: 20.00 cm  
 Height: 16.05 cm  
 Weight: 3.58 kg

**STANDARDIZED**

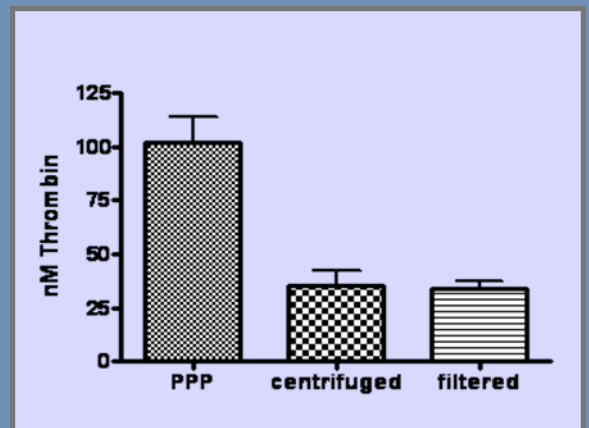

**REF 9821010 Ceveron<sup>®</sup> MFU 500**

CEV002E.002 06/2013

**Technoclone GmbH**  
 Brunner Str. 67 | 1230 Vienna | Austria  
 Tel: +43 1 86 373-0 | Fax: +43 1 86373-44  
 sales@technoclone.com | www.technoclone.com

